

Not CE Marked - Not Available for Sale in the European Union



Non-ablative
resurfacing,
redefined.

CYNOSURE[™]
LUTRONIC[™]





Achieve efficacious results with our *Integrated Cooling System*

The MOSAIC 3D integrated cooling system enables providers to safely and effectively address sensitive treatment areas around the eyes and mouth. Integrated air-cooling flows through the handpiece and onto the unique cold-compression sapphire tips. This synergy reduces pain, erythema, and edema, and is to enhance treatment efficacy, safety, and comfort.



Cold Compression Sapphire Tip



Handpiece Air-Cooling in Action






Adding new dimensions with *3D Controlled Chaos Technology™*


3D Controlled Chaos Technology (3D CCT) adds a new dimension to our patented CCT scanning.

By optimising the distribution and depth of multi-depth microbeams, 3D CCT helps reduce the risk of bulk thermal heating, allowing for an increase in treatment coagulation volume per stamp, without sacrificing safety.*

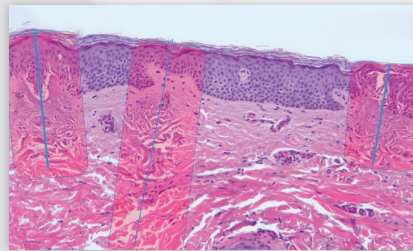

3D CCT Histology

*data on file

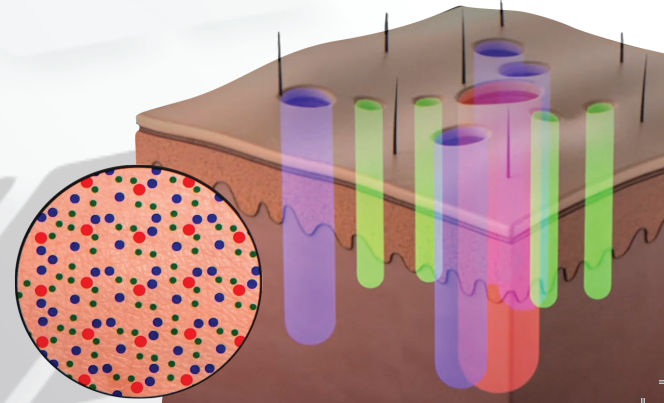


Real patients. Genuine concerns. *Life-changing results* /

Courtesy of Robert Murgia, DO



Before and after 2 treatments

Courtesy of Amy Taub, MD

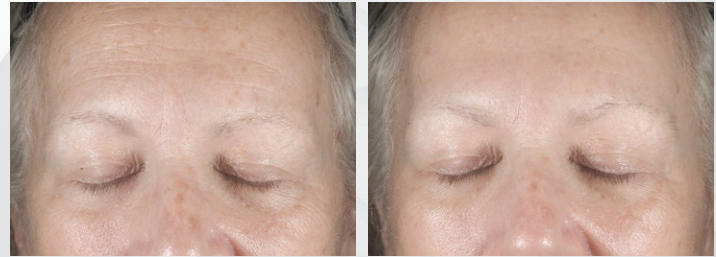

Before and after 2 treatments

Courtesy of Robert Murgia, DO



Before and after 2 treatments

Courtesy of Amy Taub, MD



Before and after 2 treatments



Distinguish your *services*

MOSAIC 3D creates a new canvas in the skin by enhancing neocollagenesis and neoelastinogenesis in the mid-dermis, leading to revitalised and brightened skin*. With its upgraded metal rollers and new sapphire tips, MOSAIC 3D rolling and stamping applications in one laser platform. Easy, air-cooled rolling enables efficient coverage and exceptional efficacy, while 3D CCT stamping in challenging areas.

MOSAIC 3D is a highly controllable and customisable fractional laser platform; its flexible treatment range enables ever-changing utility in your ever-expanding practice.

*data on file





Visit/
www.CynosureLutronicEMEA.com

Follow us:

@cynosurelutronicemea @cynosurelutronic-emea

Not CE Marked - Not Available for Sale in the European Union

Mosaic 3D is indicated for the following in CANADA: Dermatological procedures requiring the coagulation of soft tissue, as well as for skin resurfacing procedures. Indications for use include wrinkles, skin rejuvenation, acne, and various scars. Cynosure, LLC and Lutronic Aesthetic, Inc. do not diagnose or treat medical conditions, prescribe treatment protocols, provide medical advice, or otherwise practice medicine. Patient should consult with a qualified healthcare professional for diagnosis and proper treatment protocols/parameters to determine if treatment with any product is right for them. Qualified healthcare professionals are responsible for diagnosis and treatment protocols using products. MOSAIC with the 3D cube symbol represents the proprietary device name Mosaic 3D.

©2025 CYNOSURE LUTRONIC, INC. All Rights Reserved. Mosaic 3D, Controlled Chaos Technology and 3D CCT are trademarks and Lutronic is a registered trademark of Lutronic Corporation. Cynosure is a registered trademark of Cynosure, LLC.

Manufacturer of Record:
Cynosure Lutronic
5 Carlisle Road Westford, MA 01886 USA
800-886-2966



משווק בלעדית ע"י אמי טכנולוגיות
09-7760300
ami.co.il | amiweb@ami.co.il



מק"ט: CAT_MOSAIC_0426_DR_EN

Not CE Marked - Not Available for Sale in the European Union

Attributes & Specifications

Attributes	Specifications
Laser Type	1550nm, Er: Glass
Laser Delivery System	Optical Fiber with Handpiece
Pulse Energy	4 to 70mJ
Pulse Rate	up to 447.4Hz
Cooling System	Integrated Air Cooling
Cold Compression Sapphire Tips	Circle: up to 12mm diameter Tx area Rectangle: up to 8x20mm Tx area Square: up to 20x20mm Tx area
Controlled Chaos Technology™	Exclusively available in 2D and 3D
Dimensions & Weight	16" W x 23" D x 69" H & 73lbs (33kg)
Electrical Requirements	AC100-240V, 50Hz/60Hz (Fuse AC125V/6.3A) Max Power: 350VA

