

tixel[®]



Protocol



Step 1

Before treatment

Wash treatment area with soap and disinfect skin with chlorohexamin.

Step 2

Review patient's skin condition and skin type.
Refer also to clinical guidelines in user manual

Step 3

Mark face with dermatological pencil

Treatment Protocols

Step 3



Step 4

Perform treatment according to following guidelines

Skin rejuvenation - Severe wrinkleless

Skin Rejuvenation

Skin types I-III

Suggested Settings (Max.)

High Ablative

Single pulse

20% overlap

No Emla

Moisture 4 X day, for 4 days

Makeup after 2 days with sun block

Use these settings to treat cases of **severe** wrinkles and photo damage

Recommended 4 treatments, 4-5 weeks intervals

*Use periorbital HP, 6 sessions

Maintenance every 3 months



Skin rejuvenation - Mild wrinkleless

Skin Rejuvenation

Skin types I-III

Suggested Settings (Max.)

Boarder Line

Single pulse

20% overlap

No Emla

Moisture 4 X day, for 4 days

Makeup after 2 days with sun block

Use these settings to treat cases of **mild** wrinkles and photo damage

Recommended 4 treatments, 4-5 weeks intervals

*Use periorbital HP, 6 sessions

Maintenance every 3 months



Treatment Protocols

Skin Rejuvenation - Skin types V-VI

Skin Rejuvenation Skin types V-VI

Suggested Settings (Max.)

Non Ablative

Single pulse

20% overlap

No Emla

Moisture 4 X day, for 4 days

Makeup after 2 days with sun block

Recommended 4 treatments,
4-5 weeks intervals

*Use periorbital HP, 6 sessions

Maintenance every
3 months

PIH prone skin: apply HQ 4%, 3 weeks
before Tx and sunblock 3 times/day
Repeat 2 weeks post Tx



Acne Scars & Stretch Marks

Suggested Settings (Max.)

Single pulse

20% overlap


No Emla

Moisture 4 X day, for 4 days

Makeup daily after 2 days
with sun block

Recommended 3-5 treatments, 4-5
weeks intervals

Maintenance every
3 months

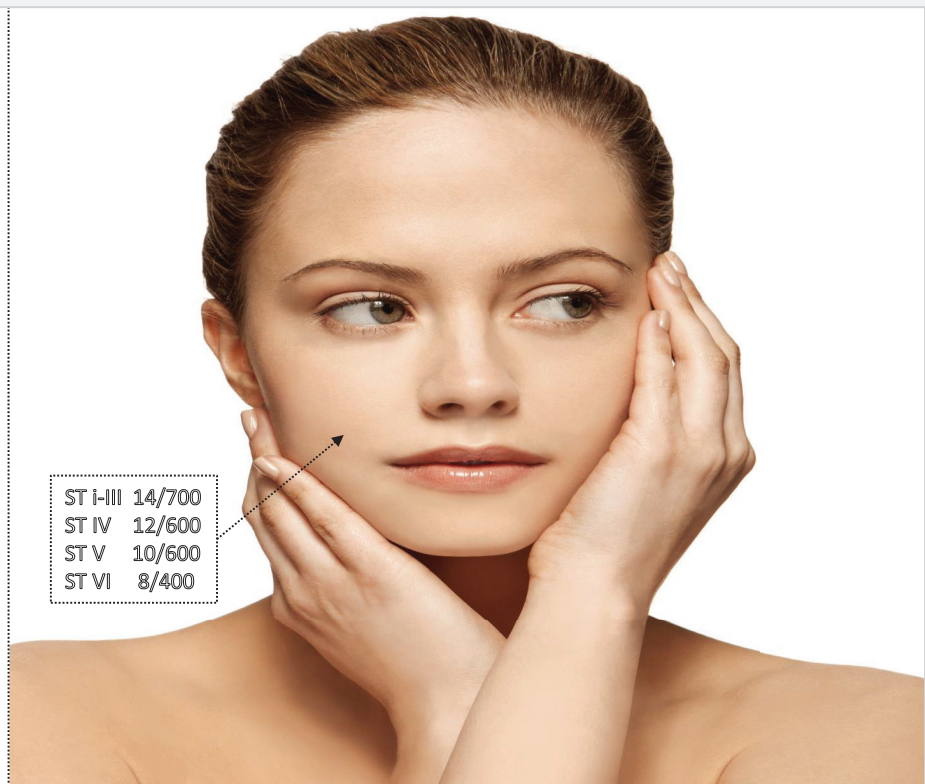

Open Channel

all skin types
20% overlap
No Emla

Approved* injectable
formulations may be applied
topically after
treatment.

Let skin rest a few hours.
Moisture 4 X day, for 4 days
makeup daily after 1 day
with sun block

Recommended 3-5 treatments,
4-5 weeks intervals

Maintenance every
4-6 months

* By local authorities



6/400
Low Non Ablative

התוויות נגד

מטופלים עם התוויות הנגד הבאות אינם יכולים לעבור את הטיפול בטיקסל:

- כל מטופל/ת שעבר/ה שיזוף במהלך 4 השבועות לפני הטיפול ו/או כל מטופל/ת המתכנן/ת לעבור שיזוף במהלך 4 השבועות לאחר הטיפול (מטופלים אשר עלולים להיחשף לשמש במשך פרקי זמן קצרים אינם מהווים קונטראינדיקציה, כל עוד הם נמרחים בקרם הגנה עם מקדם SPF50 ומעלה).
- אלרגיה לזהב.
- התפרצות של הרפס.
- סרטן העור, נגעים ממאירים ו/או נגעים טרום ממאירים מתקדמים או שומות באזור הטיפול.
- אישה בהריון או מתכננת להיכנס להריון בתקופת הטיפול, או אישה מניקה.
- אנשים עם מערכת חיסון לקויה או אנשים נוטלים תרופות מדכאות מערכת חיסונית.
- אנשים עם הפרעות ביצירת קולגן, היווצרות KELOID ו/או חוסר ריפוי נורמלי של פצעים.
- אנשים שעברו טיפולים פולשניים באזורים שמתוכננים לטיפול בטיקסל בשלושה חודשים הקודמים לטיפול עם טיקסל, או כאלה שהפצעים טרם החלימו לאחר כל טיפול אחר.
- כל אדם הנוטל או שנטל בעבר תרופות (לרבות תכשירים לשימוש חיצוני), טיפול צמחי מרפא (אוראלי או תכשירים מקומיים), תוספי מזון או ויטמינים, אשר עשוי לגרום לעור רגיש או לפגוע בתהליך ריפוי העור במהלך 3 החודשים האחרונים.
- כל אדם שהשתמש בתרופות הבאות באופן אוראלי ISOTRETINOIN® (ROACCUTAN או ACCUTANE) ב- 3 חודשים אחרונים לפני הטיפול.
- כל אדם שיש לו היסטוריה של כשל ביכולת יצירת קרישי דם או נוטל תרופות נגד קרישיות דם.
- כל אדם שיש לו קעקועים או איפור קבוע באזור המטופל.
- כל אדם שיש לו כוויות עור, יבלות בעור, עור מגורה או עור רגיש בכל האזורים שיטופל בטיקסל.
- כל אדם שעבר הרמת חוט של האזור להיות מטופלים ב- 3 החודשים האחרונים.



tixel®

NoVoxel®



Technologies

09-7760300



ami.co.il

amiweb@ami.co.il



משווק בלעדית ע"י אמי טכנולוגיות חפשו אותנו ב: